



**FRESENIUS
KABI**

caring for life



CompoLab TM

Hemoglobin analyzer

Ergonomic and easy to use, CompoLab TM is a small handheld hemoglobin analyzer with less than 2 seconds measuring time.

Fast measuring, within 2 seconds

Small size with large display

Long lasting battery capacity

Optional data management

CompoLab TM


Fast, ergonomic and easy to use

- Display always in stand-by mode – device ready to use
- No moving parts – robust for mobile use
- Light weight – approx. 180 g
- No calibration needed – minimal maintenance and cleaning
- Optional Fresenius Kabi software “HB-Master” can be installed on local PC – easy data transfer via USB cable
- Optional barcode scans – simply connect our barcode scanner to the PC
- Two charging possibilities – USB and adapter plug

Reagent free cuvette with unique design and long expiration dating

- No chemical reaction inside the cuvette – Insensitive to temperature and moisture
- Air vent prevents formation of bubbles in the cuvette
- Cuvette shelf life 2.5 years even in opened foil bag
- Can be filled from any angle
- Only one drop of blood needed

Easy cuvette
filling

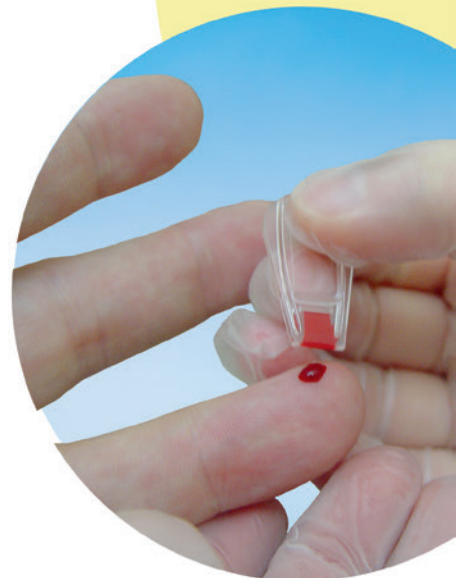

Fixed cuvette holder



Cuvette foil bag and cuvettes



CompoLab TM Hemoglobin analyzer

Process efficiency to help you achieve more



Ordering information

For more information such as literature, technical details, working instructions and equipment, please contact your local sales representative.

CE marked, Directive 98/79/EC

Electrical safety, IEC 61010-1



**FRESENIUS
KABI**

caring for life

Fresenius Kabi AG
61346 Bad Homburg
Germany

Phone: +49 61 72 608-0
Fax: +49 6172 608-5786

compolab@fresenius-kabi.com
www.fresenius-kabi.com