



A complete range of dedicated accessories

Amika Backpacks

Large


Small



Universal Table Stand



Smart Holder COM

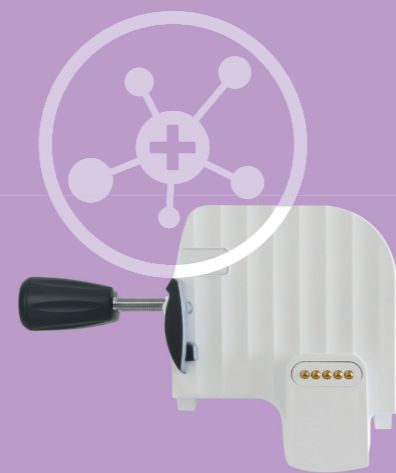

Technical details Amika+

Feeding

Continuous mode
Volume: 1 - 5000 mL (1 mL increments)
Flow rate: 1 - 600 mL/h (1 mL/h increments)

Bolus mode
Bolus Volume: 1-5000 mL (1 mL increments)
Flow rate: 1-600 mL/h (1 mL/h increments)
Bolus Number: 1 - 24, unlimited
Intervals: 15 min - 24 h (15 min increments)

Flow rate accuracy
+/- 5% at 125 mL/h

Priming
Automatic, semi-automatic and manual priming options.

Counter
Cumulative feeding volume counter from 0.001 L to 99999 L

Data event log
Feeding history: 200 events
Alarm history: 150 events

Night mode
Reduces screen and power LED brightness

Keypad lock
Disabling keypad to prevent unintentional key press

Setting lock
Option to restrict patient access to settings with access code protection.

Other settings
Activation/deactivation of target volume, flushing information, key beep, sound level adjustment; contrast & brightness settings; time between alarm sounds setting; time between target volume almost reached message and target volume reached setting

Intended use

For enteral feeding use only

Information / alarms

Pump status
Administration in progress: Droplet animation
Administration stopped: Stop symbol displayed
Administration on hold (during Bolus Interval Time): Droplet animation

Information
Start reminder, flushing information, technical information available via the menu.

Alarms
Target volume almost reached, battery almost discharged, target volume reached, door opened, incorrect set installation, downstream occlusion, upstream occlusion, empty bag / air in line, empty battery, technical error. All alarms represented by light indicators, pictograms and sounds

Occlusion detection time
At 50 mL/h: < 6 min

Empty bag / air in line detection time
At 100 mL/h: 3 min maximum

Technical specifications

Pumping mechanism
Linear peristaltic pumping system

Communication via Smart Holder COM
- Serial communication protocol for interface with Electronic Medical Record (EMR)
- Nurse call system
- Connection with computer for maintenance activities

Display
Monochrome LCD 40 x 59 mm (160 x 240 pixels)

Clamp

Versatile clamp to fix on a rail (10 - 35 mm) or a pole (8 - 40 mm)

Dimensions (H/W/D) / Weight

Pump: 138 x 128 x 48 mm / 610 g
Pump Holder COM: 132 x 118 x 46 mm (without clamp) / 420 g

Battery

4.8 V / 2.2 Ah NiMH (Nickel-Metal Hydride)
Battery life: 24 h at 125 mL/h
Battery level indication on the display
Battery charging time: 6 h

Compliance standard

IEC 60 601-1, IEC 60 601-1-2 and IEC 60601-2-24
93/42/EEC Medical Device Directive and Medical Device Regulation (EU) 2017/745

Waterproofness

Pump: IP35
Smart Holder COM: IP32

Electrical compliance

Protection against leakage current: Defibrillation-proof type CF 
Protection against electric shocks: class II

Power supply

100 - 240 VAC / 50 - 60 Hz

Power consumption

Maximum 9 W in standard operating conditions

Noise

Around 25.8 dBA at 1 m (background noise 18.3 dBA)

Housing material

ABS

Amika administration set

Free-flow protection

All Amika enteral sets feature the tube clamp, preventing free-flow when the door is open or set is not installed.

Set length

Around 2.50 m, adjustable anti-free-flow tube clamp positioning

Materials

All Amika enteral sets are DEHP-free and not made with natural rubber latex

Amika accessories

Amika backpack

For ambulatory use. Exist in large and small versions

Smart Holder Power

Pump Holder Power accessory with integrated AC/DC adapter

Smart Holder COM

Pump Holder COM accessory with integrated AC/DC adapter

USB Cable

For maintenance activities. To be connected with Smart Holder COM

Nurse Call Cable

To connect the Smart Holder COM to the Nurse Call system

Partner Software

Software for maintenance activities of Amika+ pumps

Universal table stand

An alternative pole to install pump on a table

Amika AC Adapter/charge

AC adapter power supply (plug type A)

* Only compatible with Amika Holder Power and Amika Holder COM



Amika+

Enteral Feeding Pump



Fresenius Kabi AG
Eise-Kröner-Str. 1
61352 Bad Homburg Germany
Tel: +49 (0) 6172 / 686-0

Fresenius Kabi Canada Ltd.
165 Galaxy Boulevard, Suite 100
Toronto, Ontario M9W 0C8 Canada
1-877-779-7760
www.fresenius-kabi.ca





VERSATILE

- Choice of continuous or bolus feeding mode according to clinical needs
- Customizable settings to fit user preferences
- Designed for use by healthcare professionals in acute care facilities as well as by patients at home.

EFFICIENT


- Reset parameters setting function improves nurse workflow for more effective use of time
- Alarm and feed history features facilitate user-friendly monitoring of feed intake

INTUITIVE

- The intuitive design enables easy operation for all users.
- Standard graphical symbols for universal understanding
- Complete information on the screen to follow the status of the nutrition delivery at a glance



Information light indicator

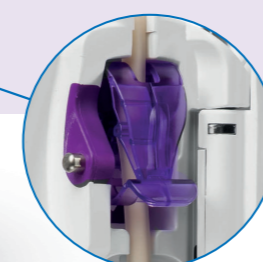


Patented anti-free flow clamp


Amika+

Personalized Nutrition Delivery

Amika+ is a user-friendly enteral nutrition pump, designed for both healthcare practitioners and patients who need advanced capabilities to effectively manage feed intake.

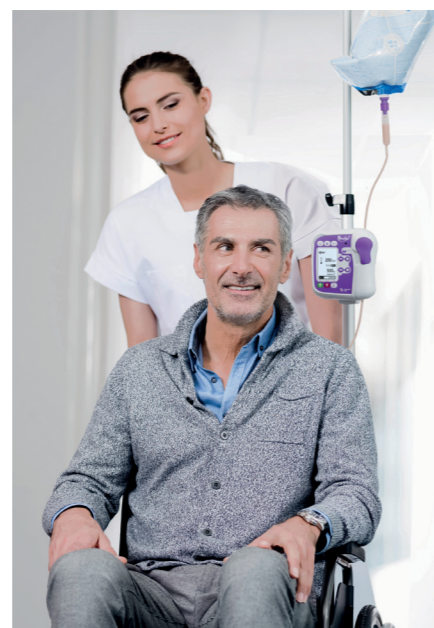


An intuitive and versatile graphical interface facilitates easy learning and use for both healthcare professionals and patients

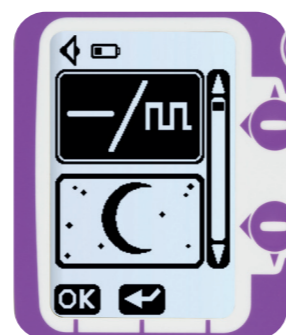

Direct access to continuous or bolus feeding mode




Reset parameters setting

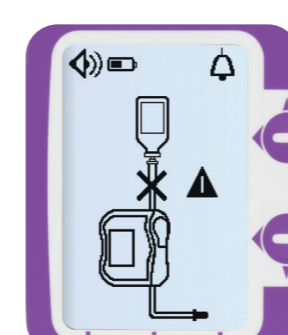

Alarm symbols for quick problem identification and action.

A Comprehensive Enteral Feeding System from Fresenius Kabi

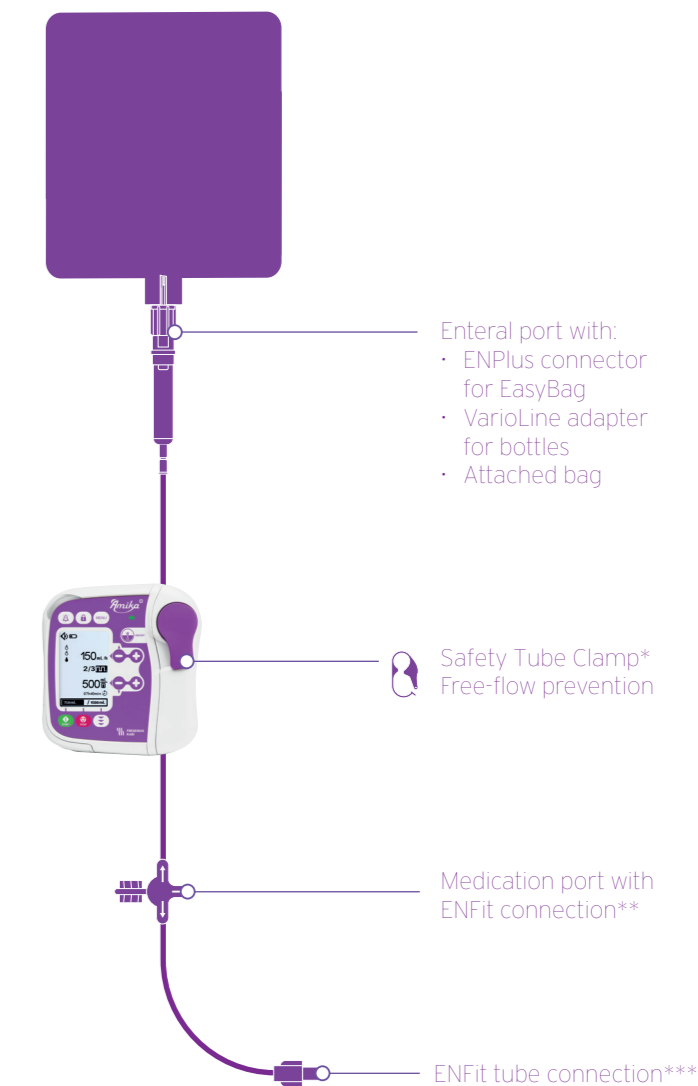

* Patented safety tube clamp

** ENFit universal sets are available with transition connectors for ENLock or LuerLock

*** ENFit universal sets are available with transition connectors for ENLock, LuerLock or funnel feeding tubes

